

FreedomAlert[®]

Enjoy the Freedom of Independent Living



INSTALLATION & OPERATION INSTRUCTIONS

Model #35911



PLEASE READ THIS MANUAL COMPLETELY AND SAVE IT FOR REFERENCE.

Thank you for choosing FreedomAlert® – the most unique Personal Emergency Response System in the world, and covered by both US and Canadian patents.

FreedomAlert is your personal link to Family, Friends, or Emergency Personnel when needed. FreedomAlert allows you to program up to 4 numbers to be reached at the push of a button.

With FreedomAlert - there are:

No Activation Costs

No Contracts

No Monthly Fees

Up to 4 user programmable contacts to call anytime. You can personalize the system with a flick of a slide switch on the Base unit to:

Call up to 4 custom contacts

Call up to 4 custom contacts & then the 911 emergency operator

Or just the 911 emergency operator.

Because you are speaking through the pendant, FreedomAlert allows you to personally communicate your emergency from anywhere in your home or yard. Emergencies can happen anywhere. And just one button activation.

Traveling? Whether you are moving across town or going to visit the grandkids across the country - just unplug your FreedomAlert from the phone jack, unplug the power adaptor and throw it in your suitcase. When you get to your destination, simply plug it into the local phone jack and power outlet. No additional programming needed. And if you call 911, they will automatically have the local 'caller ID' of that location.

Voice prompts provide easy set-up and testing.

Base unit comes with 4 rechargeable batteries to provide up to 24 hours of standby power in the event of a power outage.

Need others to share the system in your home. No problem, you can add as many additional Pendants to one Base unit as you wish.

TABLE OF CONTENTS

Overview of the FreedomAlert	2
Contents of FreedomAlert kit	4
Limitation of Liability	5
Installing the Batteries	6
Connecting the Base Unit	
Power Connection.....	8
Telephone Connections.....	8
How to Handle VoIP & DSL.....	9
Location of Base Unit.....	10
Base Unit Options & Features	
Friends List and Emergency Operator	10
Programming the Unit	
Programming Telephone Numbers	11
Telephone Number Programming Worksheet.....	12
Adding Additional Pendants to the System	14
Operating the FreedomAlert	
Operating the Pendant.....	15
Hanging Up.....	15
Lanyard, Belt Clip or Wrist-Strap.....	16
Battery Check.....	17
System Check.....	18
Emergency Dialing Using the Panic Button	19
Operating Range	19
Water Resistance.....	20
Periodic Testing / Maintenance	
Replacing Batteries	19
Out of Service Storage or Transportation	19
Troubleshooting Guide	21
Disclaimer / General Warnings & Precautions / Safety	23
Warranty	24
Replacement Parts / Additional Accessories Available	27
Product Specifications	28

Contents of the FreedomAlert Kit:


The items below are included with your FreedomAlert .



Pendant



Base Unit

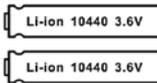

**Rechargeable Li-Ion Special
Pendant Batteries - 2 Pieces**




Rechargeable AA Batteries - 4 Pieces



Belt Clip



AC Adapter

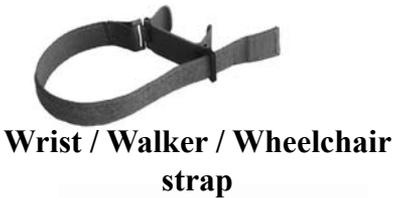

**Wrist / Walker / Wheelchair
strap**



Phone Cord



Lanyard



Screwdriver

Limitation of Liability

This page is a summary of the Disclaimer and Limited Warranty disclosed in full at the end of this instruction manual.

The purchaser agrees by using this product to the terms and conditions below and in the Disclaimer and Limited Warranty. The purchaser also agrees to read and follow all instructions and warnings on the product and contained within these Installation and Operation Instructions.

It will be the sole responsibility of the purchaser and user to assure that the FreedomAlert product is installed and programmed properly, and that the unit is used and maintained correctly. This includes, but is not limited to, periodic use to assure that the FreedomAlert, including batteries, are in proper working order, that the unit is located in an appropriate location in the home, that the electrical outlet is supplying power, and that the user has been educated as to the operation and functionality of the product as a whole.


The FreedomAlert equipment is not designed or guaranteed to prevent any loss or injury. This Limited Warranty and Disclaimer of Liability discussed in full at the end of this instruction manual constitutes the terms of sale and use of the FreedomAlert (and accessories) equipment, and if, notwithstanding these terms of sale and use of the product, there should arise any liability on the part of LogicMark LLC as a result of any cause whatsoever, regardless of whether or not such loss, damage, or personal injury was caused by or contributed to by LogicMark's negligence to any degree or failure to perform any obligation or strict products liability, such liability will be limited to an amount paid by the Purchaser for the product.

Further, LogicMark has no obligation to assure that calls are made, received or responded to, nor is LogicMark responsible for acts, or consequences of the acts, of those responding.


LogicMark provides no monitoring service for this product. It is up to the persons at the numbers called to respond in an appropriate manner.

Installing the Batteries

Pendant - small white batteries. The Pendant requires one of the white Li-Ion (Lithium Ion) batteries at any given time. The pendant should be able to run up to 2 months on a fully charged battery. The second battery is to be charging in the Base unit when not being used in the Pendant. There is no problem of over charging while in the base charger.

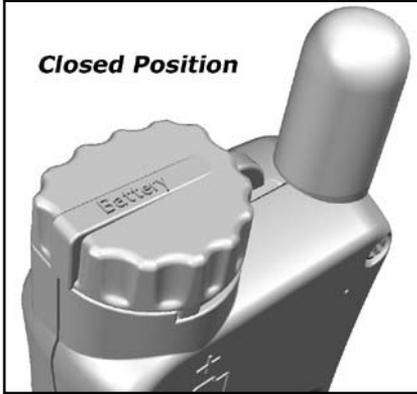

To install the battery - rotate the battery cap 90° counter-clockwise. Pull the battery cap straight up to remove. Install the battery with the positive side up facing the cap. Re-install the cap by aligning the cap as per the diagram to the left. Then push the cap back down and rotate 90° clockwise to lock in place.

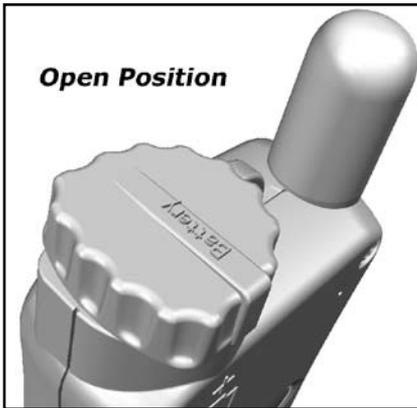


When the battery is first installed into the Pendant - you will hear an audible announcement of the battery status: "Battery is ok"; "Battery is low" or "Replace battery now". If the battery is fully discharged when you get your unit - you may not hear anything - at which point you will need to charge one of the batteries up overnight in the base charger before you will be able to use the system.


You can check the battery status at any time by simply pushing the gray battery test button on the back of the Pendant for a second.

NOTE: These white batteries are very special batteries that are not available at your local store. See Specification on page 28 for more details, or contact LogicMark to purchase replacement batteries.

Place the extra Pendant battery into the front charging area of the Base unit so that it will charge and be ready when its time to swap the batteries around. The LED on the right side will glow red when charging and green when fully charged.

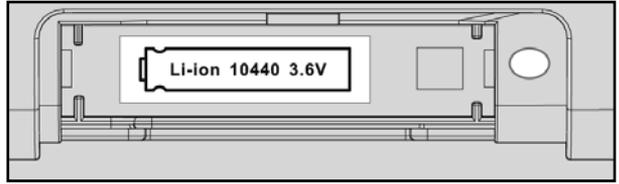

Base - larger green batteries: The Base unit requires 4 rechargeable NiMH batteries. Unlike the Pendant battery - these are readily available in many retail stores. These are to provide you with up 24 hours of battery back-up protection in the event of a power failure. They will not be damaged by continually charging in the base unit.

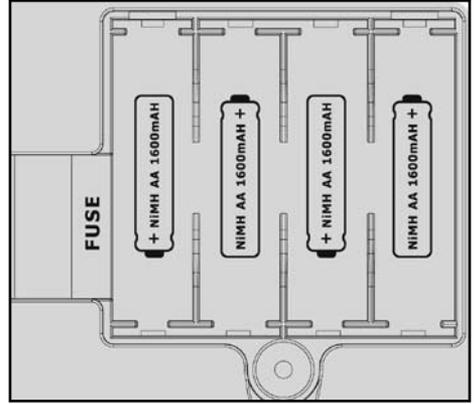


When you first get your unit - these batteries will probably be low, or even fully discharged - and will need to be charged overnight before being able to provide you the 24 hour battery back-up protection.

Unscrew the cover on the bottom of the Base unit with the supplied screwdriver. Install the 4 batteries as shown in the plastic engraving in the battery compartment area. A few seconds after you install these, you might hear “Running on battery power”. This will continue every 5 minutes for up to 15 minutes - then every 40 minutes until power is applied through the power adaptor. (If you do not hear anything - this simply means your batteries will need to be fully charged overnight. Or it may mean that you have already plugged in the power adaptor - which is why it would not announce “Running on battery power”).

NOTE: the voice prompts of “battery is OK” etc only apply to the pendant battery - not to the status of the base unit’s batteries. It is recommended that you replace the base unit batteries with a new set of similar batteries every two years.

Connecting the Base Unit

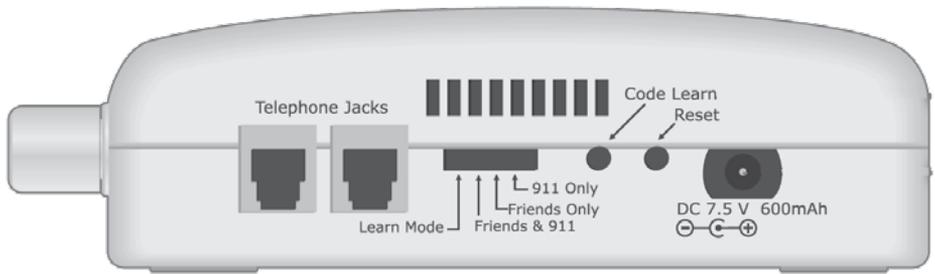

Power Connection. Simply plug in the AC Adapter in to the nearest electrical outlet. Plug the other end of the AC power adapter into the Base unit. You need to push the plug into the opening in the Base unit and twist 90° clockwise to lock in place. To remove in the future, simply rotate 90° counter-clockwise and pull outward.

Caution: Do not plug the AC Adapter into an outlet which is controlled by a switch. The switch could accidentally be turned off, thus rendering the Base unit inoperable after approximately 24 hour when the back up batteries become discharged.

Telephone Connection. No special phone line or special service is required to use your FreedomAlert. Your phone will work with your existing touch-tone phone line and existing service. It works the same way as any cordless phone. (**Touch Tone service is required**). Simply plug one end of the telephone cord into your telephone jack in the wall, and the other end into either one of the telephone jacks at the back of the Base unit. Note - you will feel a 'click' when the jacks are firmly seated in the wall and the Base unit jack recepticals.

If you need to plug a conventional telephone into same outlet used by the FreedomAlert, simply plug the telephone into the spare jack at the back of the Base unit.

How to Handle VoIP & DSL

If Your Phone System is VoIP: Voice-over-Internet Protocol (VoIP) is a telephone service that sends your call over the Internet instead of through the regular telephone system. FreedomAlert will work well on most VoIP systems.


It is important (if you haven't already done so) that you register your physical address with your VoIP provider and check to insure that 911 service is available to you.

If you are unsure if your phone system uses VoIP, or to find if your VoIP system is capable of calling 911, please contact your telephone service provider.

To install FreedomAlert, simply plug the product into the telephone jack (or adaptor jack) that you would normally plug a regular telephone into.. Note, if you have both a working traditional phone line AND a VoIP phone system we recommend connecting the FreedomAlert system to the traditional phone line.

For more information on using the FreedomAlert product with your VoIP system, visit our website at www.LogicMark.com and view the FAQ under the OUR PRODUCTS tab.

How to Handle DSL: If your home has DSL service for your computer, you will need to use a DSL filter in line with the FreedomAlert OR IT WILL NOT WORK. This filter can be obtained from a local electronics supplier or your DSL service provider. See illustration to the right.

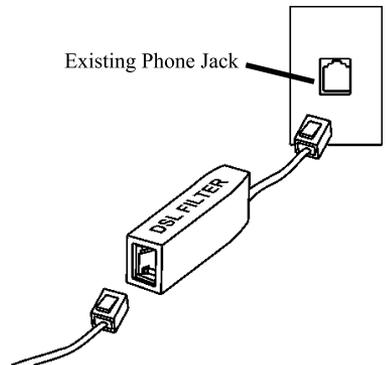

To FreedomAlert


Locating the Base Unit: For best results and longest range, locate the Base unit in an area away from household appliances which could present electrical interference. These include microwave ovens, televisions, electric mixers or blenders, hairdryers, wifi routers or other cordless phones. Do not place on a metal surface.

If more than one FreedomAlert is to be used in one home, it is recommended that the Base units be located at least 10' apart.

Base Unit Options and Features:

The FreedomAlert can work right out of the box with no programming. To avoid programming any telephone numbers, simply slide the slide switch at the back of the Base unit to the **“911 Only”** mode. Then, when you push the blue button on the Pendant - you will be immediately connected to 911 emergency services.

If you choose not to program any custom numbers, simply skip these next steps on Programming Telephone Numbers and jump to page 14.

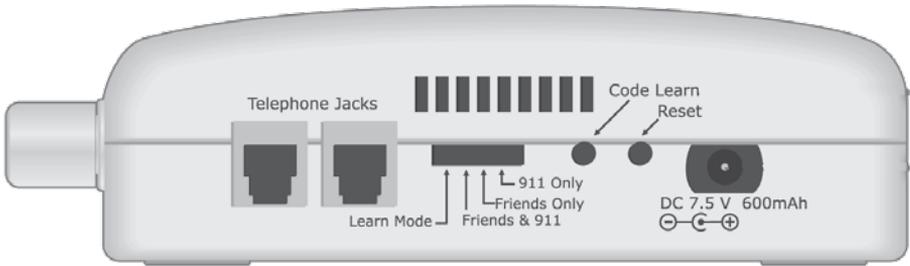

Friend's List and Emergency Operator.

In addition to being able to call 911, you can also program up to 4 telephone numbers of your choosing. You need not program all 4 - just one if you like. Typically, you might have it programmed for 2 numbers - to call your neighbor and if they are not available, to call a family member. Even if you do program a number or more into the Base unit - you can still at any time choose to call **“911 Only”** by just sliding the switch at the back of the Base unit. Your numbers remain in memory and will be ready for use when you slide the slide switch back to one of the **“Friends”** position. Note: Disconnecting power or changing batteries does not erase the numbers.

How does the Base unit know when to go to the next number? The base unit will first dial the first number in the “**Friends**” list. About 20 seconds into the conversation, the Base unit will interrupt the call and announce “*This is an emergency call. If you want to continue this call - Press 5 now. Otherwise, we will hang up and call the next contact or emergency operator*”. At this point - the person that is being called must press the “5” key on their telephone if they want to continue the call. If they don’t press “5”, then the Base unit will assume they do not want to talk, or that it may be an answering machine or a child - at which point it hangs up and goes on to the next number. It is important to instruct the people you have programmed into your “**Friends**” list to press “5” when they get this call, and when they hear this message.

Programming the Base Unit

Programming Telephone Numbers:

If you choose to only contact 911 - and not a custom telephone number - skip these next steps (and go to page 14) - and make sure your slide switch on the back of the Base unit is set to “**911 only**”.

NOTES:

- As you program the system using your telephone on the same line as you have established the call with in the previous steps, you will hear the base unit speak each number that is pushed as you push them.
- Dial slowly and listen for each number to be announced. If you did not hear the announcement - the number was not recorded.
- The FreedomAlert can learn telephone numbers up to 16 digits long, and you do not have to program all 4 numbers
- Need to add a pause in your number sequence.....simply insert a “*****” where ever you need a 1 second pause. Insert “******” if you need a 2 second pause.
- Remember to program a **1** before any long distance numbers.
- Remember to program any **AREA CODE** numbers if needed.
- Remember to program any **PREFIX NUMBERS** (like 9 to get an outside line) if needed.

Step 1: Write out the numbers you are wanting to program into the system. Below we will refer to Number 1 as the first number, Number 2 as the second number etc. Include area-codes and a 1 or 9 if they are needed for dialing.

Step 2: Establish a phone connection with a friend or family member - or with your cell phone. This is done as follows. From your home phone (on the same line as the FreedomAlert is hooked-up to) call your cell phone (or a friend or neighbor if you don't have a cellphone). Answer the phone (or wait until your neighbor answers the phone and tell them just to wait a minute while you program your system). Note: If you don't do this step, then as you begin to program the unit with your telephone - the phone company will think you are trying to dial an outline number and will try to complete the call. By having the phone line already in use with your cell phone or friend or family member - you are now ready to program the FreedomAlert.

Step 3: With the phone connection established in Step 2, slide the slide switch at the back of the Base unit to “**Learn Mode**”. You will hear an audible announcement of this position.

Step 4: Determine how many contact numbers you want to enter. Follow the next steps to enter these numbers.

NOTE: The “#” key (the key in the lower right corner of your phone keypad) must be entered before each phone number as shown in the work sheet on the next page - and then pressed 4 times to exit out of programming.

Telephone Number Programming Worksheet:

Note: You will be using your telephone keypad on an open telephone line to program the emergency call telephone numbers. The base will speak each number that is pushed on the telephones keypad.

Dial Slowly – listen for each number to be announced, if you did not hear the voice announcement – the number was not recorded.

The Freedom Alert can learn telephone numbers up to 16 digits long and you DO NOT have to program all 4 contact telephone numbers. When you have finished entering the telephone numbers that you want, skip to step 9.

1. Completely fill out this worksheet before starting to program telephone numbers.
2. Call a friend or family member on your telephone.
3. Slide switch to LEARN MODE [voice prompt – Learn Mode]
4. Enter: 1234 [voice prompt – one, two, three, four, programming mode]
5. (Enter “#” then the 1st phone number) :

[voice prompt – numbers as you enter them]
Skip to step 9 to Exit Programming, or
6. (Enter “#” then the 2nd phone number)

[voice prompt – numbers as you enter them]
Skip to step 9 to Exit Programming, or
7. (Enter “#” then the 3rd phone number) :

[voice prompt – numbers as you enter them]
Skip to step 9 to Exit Programming, or
8. (Enter “#” then the 4th phone number) :

[voice prompt – numbers as you enter them]
9. Enter: ##### [voice prompt – end of programming]
10. Slide switch back to Friends & 911, Friends Only or 911 Only
11. Programming complete!

Adding Additional Pendants to the System.

FreedomAlert comes with one Pendant. If you want to add additional pendants to your system - follow the directions here. If not, skip this section.

The Pendant's unique ID code should already be learned into the Base unit when it came from the factory. You can add additional Pendants (up to 16) to the system by simply teaching each pendant's ID code to the Base unit.


To program a new Pendant to the Base unit:

1. Push the RED Code Learn button on the back of the base unit. You will hear "Synchronizing Mode"
 2. Immediately PUSH and HOLD both the Blue panic and Gray battery test buttons at the same time. You will hear "Synchronizing" from the pendant.
 3. Release both button on the pendant.
 4. The pendant is now learned to the base.
 5. If teaching additional pendants, repeat steps 1-3
- The base unit will automatically exit the pendant learning mode after several seconds.
 - Teaching additional pendants does not erase previously learned pendants.
 - If you hear the voice prompt "Base and Pendant Out of Range" an error has occurred, wait 1 minute for the base to reset and repeat steps 1-4.

At a later point in time, you might want to erase your previously learned numbers. You can erase all previously learned numbers by simply holding down the red **Code Learn** button at the back of the Base unit for 8 seconds. you will hear the confirmation of "Previous Pendant Codes Erased" from the Base unit. After that, you can program a new pendant to the base unit following the steps above as desired.

Operating the FreedomAlert:

Operating the Pendant. The Pendant is simple in operation - having one main Blue button - and a smaller Gray button on the back.

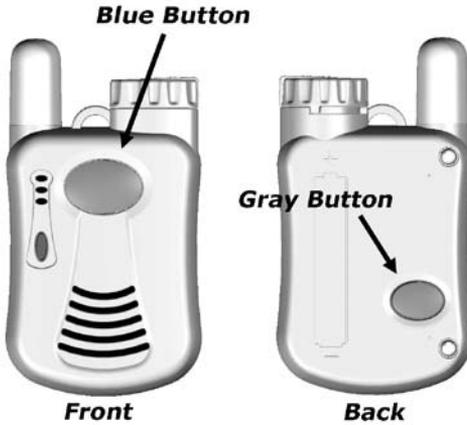

The Blue button on the front is the main button to push for a full second to initiate a call. Once this button is pushed the call proceeds forward to the 'Friends' list and or 911 depending on the slide switch setting on the back of the Base unit.

Holding the Pendant. At this point - you simply hold the Pendant as you would any normal telephone receiver - close to the head so that you can hear the other person being called - and they can hear you speaking.




Hanging -Up. If you choose to cancel the call once initiated - you can simply push the Gray button on the back for a second to cancel the call. you will hear "Hanging up". If at the end of the conversation, the other side hangs-up, the system will hang-up automatically over the next minute or two. Or, you can push the Gray button to end the call immediately.

Lanyard, Belt Clip or Wrist Strap.

FreedomAlert comes with 3 accessories for carrying the pendant . You can carry the Pendant using the Lanyard. You can slide the Pendant in the Belt Clip, or you can use the Wrist-strap as a means to carry the Pendant. Note that the Wrist-Strap can also be used as a way to attach the Pendant to a wheel chair. Choose the method that suits you best.

The drawing at the right shows the method of attaching the lanyard to the Pendant.

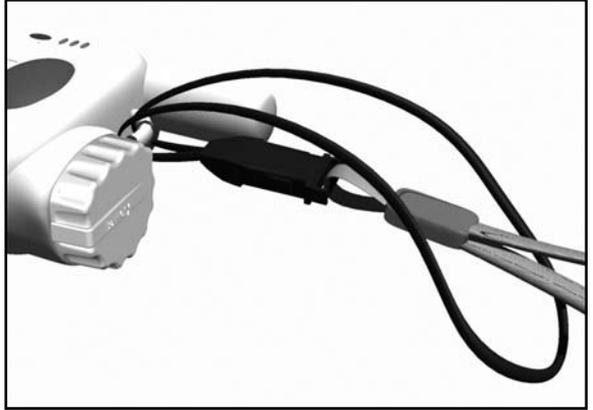

Battery Check:

The FreedomAlert has a built-in battery tester. When you press the Gray button at the back of the Pendant, a voice announcement will inform you of the battery condition. One of the following messages will be heard:

“Battery is ok.”

Check the battery condition twice a month.

“Replace battery soon.”

The battery should be replaced within the next week.

“Replace battery now.”

The battery should be replaced now.

To replace the battery - simply exchange the battery in the Pendant with the battery that has been charging in the Base unit.

Please note: If there is no response after pressing the button for one second or more, the battery is dead and must be replaced immediately.

The battery should be tested at least twice a month along with the System Check outlined on the next page.

When it gets to the point after a few years that the pendant battery only lasts 3 or 4 weeks before needing to be charged - it is then time to replace both the Pendant batteries. Note: This is a very unusual battery not available at most outlets. Order replacement batteries from LogicMark as directed on page 27.

System Check:

The button on the bottom of the pendant case which is used to perform the BATTERY CHECK, will also perform a SYSTEM CHECK. Press and hold this button for more than four (4) seconds. One of the following announcements will be made.

“System check, all systems are ok.”

This confirms the following:

Battery status.

The Pendant is working.

The Base unit is working.

The Base unit is connected to a working phone line.

The Pendant is within range of the Base unit.

“System cannot detect dial tone.”

The Base unit is not connected to the phone line or can not get a dial tone for some reason.

“Base and Pendant out of range.”

This confirms one or more of the following:

The Pendant has malfunctioned.

The Base has malfunctioned.

The Pendant and Base unit are not communicating with each other or are out of range.

If you receive a failure notification, please refer to the TROUBLESHOOTING guide on page 21 of this manual to determine the cause.

This testing feature is very useful - allowing you to confirm that the system has coverage from all parts of your home. Simply walk to the various parts of your home and push the Battery Check / System Check button until you hear the dial tone. You can then release the button, and you will hear “System Check. All Functions are OK” voice prompt, or “Pendant and Base out of Range”.

We highly recommend that you perform this System Check at least twice a month to confirm that your system is working properly and that the pendant and base unit are in good communication and that the base unit can get a dial tone. It is also recommended to perform this system check after a lightening storm that could damage the base unit which is plugged into the phone line.

Emergency Dialing Using the Pendant:

In the event of an emergency, press the blue panic button on the front of the Pendant. When this button is pressed for more than 1 second, the FreedomAlert will automatically start dialing from the list of numbers.

If set to the “911 only” mode, then the 911 emergency operator will be called. When the emergency response operator answers, speak in a normal voice and you will be heard.

If set to the “Friends only” mode - the system will dial the first telephone number in the list. If no one is available for that number, it will dial the second number, etc - until it gets to the end of the list. It will then repeat from the top of the list.

If set to the “Friends and 911” mode - the system will dial the first telephone number in the list. If no one is available for that number, it will dial the second number, etc - until it gets to the end of the list. After that, it will then dial the 911 emergency operator.

Operating Range:

The operating range of your FreedomAlert will depend on many factors including intervening walls, electrical interference or various appliances such as vacuum cleaners, microwave ovens, mixers, coffee grinders, hair dryers and other sources of electrical noise around the house. The system should cover your typical home and a short distance into the front, side and back yards. Perform the System Check on page 17 to determine the boundaries of your system to know the limits of operation in your home and surroundings.

Water Resistance:

FreedomAlert will resist occasional splashes from water - i.e as in the shower. However, it should never be submerged or exposed to water continuously.

Periodic Testing / Maintenance

It is highly recommended that you preform a full System Check (see page 18) at least twice a month. This will confirm that the Pendant battery is in good condition and that the Base unit is able to get a dial tone.

Replacing Batteries: Replace the Base unit batteries every 2 years - or as needed. Be sure to purchase NiMH AA batteries with 1,500 mAh or greater capacity. A set of batteries in good condition and fully charged will supply 24 hours of back-up in the event of a power outage. You can test this by unplugging the AC Adapter and hearing the voice prompts. The Base unit will announce “Running on battery power” once every 5 minutes for one hour - then once every 40 minutes until the batteries are depleted.

Out of Service Storage or Transportation: When storing the Base unit -or taking it on a trip - or sending it in for repair - it is recommended that you remove one or all of the AA rechargeable batteries from the base unit so that the “Running on battery power” announcement for the next 24 hours doesn’t drive the postmaster crazy. Remember to re-install the 4 batteries in the Base unit when it is put back in service.

TROUBLESHOOTING GUIDE

What if I can't get the FreedomAlert to work?

Step 1 – Push the gray RESET button with a pen at the back of the base unit (next to the power cord). Check the system. If it is now working skip the remaining steps.

Step 2 – Check the PENDANT.

What to do:

Press the Gray Battery Test / System Check button on the back of the Pendant.

What to expect:

A voice should announce the condition of the battery.

A light on the front of the Pendant will come on.

If nothing happens, the most likely cause is that the battery is dead or has been installed backwards.

If the Pendant says “Battery ok”, go on to step 3.

Step 3 – Check the connections on the BASE UNIT

What to do:

Make sure the AC Adapter is plugged into a live AC outlet and into the Base unit.

What to expect:

The Red POWER light will be on if the Base unit is connected to power.

If the POWER light is on and the Base unit is connected to a working phone line, go on to step 3.

Step 4 – Make sure the Pendant and the Base unit are communicating with each other.

Background information: The Pendant and the Base unit are connected by

a wireless radio link. They must be synchronized to work together. They should have already been synchronized at the factory, but resetting this is simple.

What to do:

Push the red Code Learn button at the back of the Base unit for a second. You will hear “Synchronizing Mode”. Then within a few seconds - PUSH and HOLD both of the Blue panic button and the Gray battery test button on the back TOGETHER for a few seconds. You will hear “Synchronizing” from the Pendant. Release both buttons. You now have your Pendant synched to your Base unit. You can add as many new Pendants as you like to the system at any time (up to 16) - without having to worry about the previous ones being erased.

Now perform the System Check test again (page 18) to make sure all is well.

If your FreedomAlert is still not working, call our TECH SUPPORT department at 703-934-7934 (option 2 for customer service).

Information The FCC Wants You To Know:

FCC ID: TYD35911

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Privacy of communications may not be ensured when using this product.

Canada IC: 8471A-35911

General Warnings And Precautions!

Warning!

- LogicMark, LLC DOES NOT represent this unit to be waterproof. To reduce the risk of damage to the unit, DO NOT expose the pendant or base unit to prolonged water, rain or moisture.
- FreedomAlert is intended for residential use only. It may not work properly if connected to some commercial phone systems.
- FreedomAlert works with tone dialing systems only. Do not use if you have pulse dialing only.
- This product and/or accessories may contain lead, a chemical known to the State of California to cause birth defects or other reproductive harm.

Note: FreedomAlert does not interfere with pacemakers.

Important Safety Instructions:

When using your FreedomAlert , please observe the following instructions in order to reduce the risk of personal injury, electrical shock, or fire.

1. Save these instructions for future reference.
2. Read all instructions carefully and make sure you understand them.
3. Unplug this product from the electrical outlet before attempting to clean it. Do not use any liquids for cleaning either the Base unit or the Pendant. Use only a soft damp cloth.
4. Do not locate the Base unit near any source of water, such as in the bathroom, near a sink in the kitchen, or the like.
5. Do not open the unit, (except the battery compartments as directed in these instructions) on the Pendant for any reason. Any attempt to open the case of the Base unit or the Pendant will void the warranty and may damage the unit or prevent it from working properly.
6. Do not use any power source other than that supplied from LogicMark, LLC specifically for this product. Use of any other power source could damage your unit or make it unsafe.
7. Always dispose of batteries properly. Never thrown batteries into a fire. Consult your local ordinances for proper battery safety.
8. If your unit becomes damaged do not attempt to repair it. Please return it to LogicMark, LLC for service or replacement as stipulated on the next page under “ONE YEAR LIMITED WARRANTY”.

LogicMark Disclaimer and Limited Warranty

DISCLAIMER: LogicMark, LLC. will not be responsible for the improper use of FreedomAlert, nor will it be responsible for failure resulting from the use of other equipment connected to the same phone line. We will not be responsible for the quality of the phone line or the reliability or quality of the phone service with which FreedomAlert is used. LogicMark, LLC. will not be responsible for the installation of the FreedomAlert. It will not be responsible for the improper use or abuse of this unit.

COMMUNICATION AND RESPONSE LIMITATIONS: Purchaser acknowledges that signals which are transmitted over telephone lines, or other modes of communication pass through communication networks wholly beyond the control of LogicMark LLC and are not maintained by LogicMark LLC, and, therefore, LogicMark LLC shall not be responsible for any equipment or communication failure which prevents transmission signals from reaching your contact list including emergency 911 operators or damages arising therefrom. Purchaser acknowledges that LogicMark LLC provides no response to its FreedomAlert equipment. The equipment is designed to communicate with the contact list of your choice as well as emergency 911 operators and LogicMark LLC is not and shall not be responsible for ambulance, police or other emergency response time or that any response will be provided by the anyone on your contact list including the emergency 911 response center.

ELECTRIC AND TELEPHONE SERVICE IS NECESSARY AND PURCHASER'S RESPONSIBILITY: Purchaser acknowledges that the FreedomAlert equipment (including optional Battery Back-up accessory) plugs into a standard telephone jack and communicates over standard telephone lines using two way voice communication. You must use a special adaptor (not supplied) for the equipment to work with VOIP Internet transmission. VOIP Internet will not work when your electricity service is interrupted unless you have battery back up for your computer. Purchaser acknowledges that the FreedomAlert equipment requires 110 Volt AC power and electrical outlets and receptacles, telephone hook-ups, RJ11 Block or equivalent, all of which is Purchaser's responsibility to obtain and maintain.

TESTING AND SERVICE OF FreedomAlert EQUIPMENT: The FreedomAlert equipment, once installed, are in the exclusive possession and control of the Purchaser, and it is Purchaser's sole responsibility to test the operation of equipment and request warranty service if the equipment is under warranty.

LIMITED WARRANTIES OR REPRESENTATIONS. LogicMark LLC does not represent nor warrant that the FreedomAlert system will prevent any loss, damage or injury to person or property, or that the personal emergency alert system will in all cases provide the protection for which it is installed or intended. Purchaser acknowledges that LogicMark LLC is not an insurer, and that Purchaser assumes all risk for loss or injury to Purchaser's property or person. LogicMark LLC has made no representation or warranties, except those expressed herein and hereby disclaims any warranty of merchantability or fitness for any particular use.

PURCHASER'S EXCLUSIVE REMEDY: Purchaser's exclusive remedy for LogicMark LLC's default hereunder is to require LogicMark LLC to repair or replace, at LogicMark LLC's option, any equipment or part of the personal emergency alert system which is non-operational during LogicMark LLC's warranty period.

LIMITATION OF LIABILITY: The FreedomAlert equipment is not designed or guaranteed to prevent any loss or injury. This Limited Warranty and Disclaimer of Liability constitutes the terms of sale and use of the FreedomAlert (and accessory) equipment, and if, notwithstanding these terms of sale and use of the product, there should arise any liability on the part of LogicMark LLC as a result of any cause whatsoever, regardless of whether or not such loss, damage, or personal injury was caused by or contributed to by LogicMark LLC's negligence to any degree or failure to perform any obligation or strict products liability, such liability will be limited to an amount paid by the Purchaser to LogicMark LLC for the product, or to the sum of \$300.00, whichever is greater.

ONE YEAR LIMITED WARRANTY

LogicMark, LLC warrants to the original consumer/purchaser that this product shall be free of defects in material and workmanship under normal use and circumstances for a period of **one (1)** year from the date of original purchase for use. When the original consumer/purchaser returns the product pre-paid to LogicMark, LLC; 8625 Hampton Way; Fairfax Station, Virginia 22039, USA within the warranty period, and if the product is defective, LogicMark, LLC will at its option repair or replace such.

This warranty shall constitute the sole liability of LogicMark, LLC concerning the product. LogicMark, LLC expressly disclaims all other warranties INCLUDING, WITHOUT LIMITATION, THE WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. NO PERSON, FIRM, OR CORPORATION IS AUTHORIZED TO ASSUME FOR LogicMark, LLC. ANY OTHER LIABILITY IN CONNECTION WITH THE SALE AND USE OF THE PRODUCT. LogicMark, LLC and agents and distributors will bear no liability whatsoever for incidental or consequential damages or charges of any kind.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above disclaimer regarding incidental or consequential damages may not apply to you.


This warranty is void if the product has been damaged or tampered with or if the product or any such parts have been opened. In all cases of damage during shipment, a claim must be filed with the shipping carrier and not with LogicMark, LLC.

This warranty gives you specific legal rights; you may also have other rights which vary from state to state.

OUT OF WARRANTY REPAIRS

If the warranty period has expired or if you are not the original owner of the product, LogicMark, LLC will at its option either (1) replace this product with a functionally similar (but not necessarily identical) refurbished product or (2) repair the original product and return it to the original consumer/purchaser after payment of repair/replacement charges has been received.

Replacement Parts / Optional Accessories:

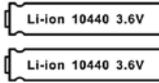

**Rechargeable Li-Ion Special
Pendant Batteries 2 Pieces**
Part #35917 \$15.00



**Rechargeable AA Batteries
4 Pieces**
Part #35918 \$15.00



Lock Box
Part #30913 \$34.95



Extra Pendants
Part #35915 \$119.95

All prices include shipping and handling to the 50 United States. Add an additional \$7.00 for Canadian shipments. Orders can be sent in with a check or credit card information. Include your full mailing address and telephone number. Orders can also be placed via telephone at 1-800-519-2419.

Always include your pendant and base serial # (on the bottom of the case) when ordering. Or have it ready when you call in.

Product Specifications:

RF Characteristics:	900 Mhz Spread Spectrum
Operating Range	Covers your typical American house and into the front, back, and side yards.
Power Adapter Ratings	Input: 110 - 240 VAC Output: 7.5 volts DC - 600 mA
Base Unit Power Consumption	50 mA in standby mode 200 mA when dialing
Back-up Battery Supply	AA NiMH 1,500 mAh batteries (4 pieces). Should be replaced every 2 years.
Back-up Operation Duration	24 hours with fresh batteries
Base Unit Fuse	5 x20 mm fuse. 250v, 0.5 amp (Radio Shack Part #270-1061)
Dialing Style	Touch-Tone only
Telephone Number Length	16 digits maximum
Telephone Numbers	Can learn up to 4 numbers plus the emergency number.
Pedant water-resistance	Water Spray Standard IPX4 (shower only).
Number of Pendants/Base Unit	16 Pendants can be programmed to one Base unit.
Pendant Batteries	3.6 (or 3.7) volt 500 mAh Li-Ion size 10440. (2 pieces) Should be replaced very 2 years Purchase replacements from LogicMark
Operating Temperature	32° to 120° F. (0° to 49° C).

LogicMark

8625 Hampton Way;
Fairfax Station, Virginia 22039

Tel: 1-703-934-7934

Toll Free: 1-800-519-2419

Fax: 1-703-934-7935

www.LogicMark.com

FreedomAlert is covered by the following patents:
US # 5,521,582 & Canadian # 2,080,921
Other patents pending

©2009 LogicMark, LLC